



Ease of Wound Care at Home

TOPIVAC TOPI

MOBILE DISPOSABLE NPWT

WOUND THERAPY DEVICES & SETS

**5 CLASSIC DRESSING SET INCLUDING
TOPI A100 MOBILE DISPOSABLE DEVICE**



3 different vacuum modes;

**-75
mmHg**

**-100
mmHg**

**-125
mmHg**

It's a negative pressure wound treatment with a single treatment protocol, designed for single patient use. The negative pressure wound closure system accelerates tissue healing by increasing microcapillary circulation. The device and dressing set operate in harmony and integration with each other.

The system is an active wound treatment system that maximizes patient treatment comfort. Its compact and ergonomic mobile structure (the size of a mobile phone) offers patients quality life and comfort of movement without restricting their daily lives and mobility during wound treatment processes.

While it heals wounds quickly and effectively, reduces the need for expert personnel with its ease of application and provides ease of home wound care and mobility to patients. Reduces hospitalizations, lowers treatment costs and provides more efficient and widespread delivery of health services.

**5 INCISION DRESSING SET INCLUDING
TOPI A200 MOBILE DISPOSABLE DEVICE**



3 different vacuum modes;

**-35
mmHg**

**-50
mmHg**

**-65
mmHg**


TOPI A100 and A200 continuously draw exudate from the wound via a special dressing set applied to the wounded area specially produced for the device; provide equal distribution of negative pressure to the wounded area. Heal the wound by forming granulation. Unlike the A100, the A200 is designed for use post-operative incision wounds. Vacuum values are lower and with its special medical foam specific for incision wounds, keeps the incision edges together without damaging the sutures, provides rapid healing of the wound and helps to prevent scarring aesthetically.



Diabetic Foot Wound



Chronic Wound

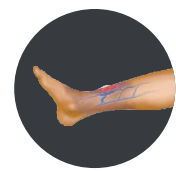

Venous Ulcer



Pressure / Bedsore Wound



Burn Wound



Infected Wounds

TOPI A100 CLASSIC


Used for rapid healing of all kinds of wounds up to 100 cm² with the approval of the physician.

WHICH WOUND TREATMENTS IS IT USED?

TOPI A200 INCISION



Plastic Surgery Wounds

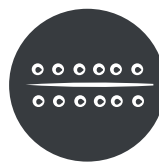

Caesarean Section Wounds



Incision Wounds



**Cardiothoracic & Vascular
Surgical Wounds**



**Obstetric & Gynaecological
Surgery Wounds**



**Orthopaedic
Surgery Wounds**

Used by physicians after surgery to reduce the risk of hospital infection, to prevent the wound from becoming infected, to heal the wound quickly and to prevent scarring.

TOPI Mobile Disposable NPWT Vacuum Assisted Wound Care Therapy Set


TOPI Mobile Disposable NPWT Vacuum Assisted Wound Care Therapy Set

Includes 1 Protocol and 5 Dressing Treatment Sterile Products.

- 1 TOPI Mobile Disposable NPWT Vacuum Assisted Wound Care Therapy Device (TOPI A100 or TOPI A200)
- 5 Pcs Dressing Foam Set (Classic or Incision)
- 5 Pcs 75 cc Canister
- 1 Pcs Extension Hose
- 1 Pcs Extra 30x40 cm Drape
- 1 Pcs Device Carry Bag
- 1 Pcs User Manuel



TOPI A100 Classic NPWT Vacuum Assisted Wound Care Therapy Device

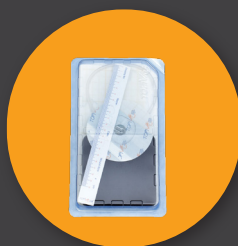

Classic Dressing Foam Set
10*10 cm Dressing Foam
1 Pcs Trackpad
1 Pcs 30x40 cm Drape
1 Pcs Ruler



30*40 cm Drape



Extension Hose



TOPI A200 Incision NPWT Vacuum Assisted Wound Care Therapy Device



Incision Dressing Foam Set
4,5*22 cm Dressing Foam
1 Pcs Trackpad
1 Pcs 30x40 cm Drape
1 Pcs Ruler



75 cc Canister



Device Carry Bag

ADVANTAGES



Practical and Easy to Use

- Ease of starting and ending treatment with one button.



Visual & Audible Warning System

- High security with visual & audible warning system.
(Instant alert during vacuum leak, full canister and low battery)



User Friendly

- Zero cross-contamination risk.
- Zero cost for preventive maintenance.
- Canister-based system for effective exudate management.
- A small lightweight design that is easy to wear or hang on the body, can be hidden under clothes, **helps you continue your social life uninterrupted without anyone noticing.**
- Thanks to its innovative design that allows silent operation during therapy, does not disturb you even in your sleep.
- Uninterrupted treatment thanks to efficient battery usage.



Ease of Use Without Expert Personnel

- Adaptable from hospital to long-term care, outpatient settings or home care.



Skin Friendly

- Optimum breathability for the patient's healthy skin while maintaining a vacuum environment for the wound.
- Special foam system used for sucking exudates and any wound fluids by creating a hygiene environment in the wound.



Movement Comfort

- Special dressing set provides a strong adhesive bond to maintain an active lifestyle while keeping the dressing intact during treatment.
- Quality life and comfort of movement with its compact mobile structure.

Beyond That Is The Patient Comfort...

Advanced Wound Management System

TOPIVA[®] TOPICAL
WOUND CARE VACUUM
TREATMENT
UNIT



topi.vac



Topivac NPWT



TopiVac Advanced Wound Care



TopiVac - @topivacwoundcare



+90 (312) 354 60 40



info@amedus.com



www.topivac.com



Ostim OSB Mah. 100. Yıl Blv No: 99/34

Yenimahalle 06374 Ankara/ Turkey

